



Competence CP323S

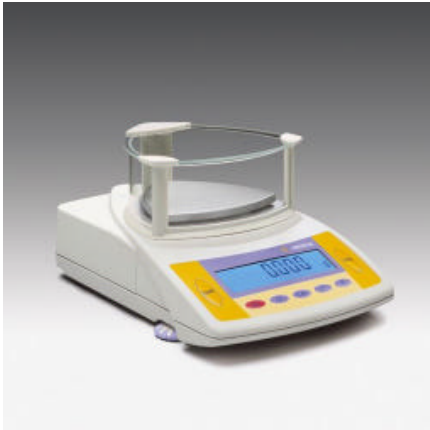

Key features. The optimum for all weighing tasks.

- Practical draft shield: hinged cover for easy sample loading from the top; removable side panels for precise filling; ensures thorough cleaning
- Excellent, easy-to-read backlit display with especially large digits (16.5 mm)
- Monolithic weigh cell technology – rugged, precise, reliable – year after year
- Very easy to operate; all keys feature positive click action; 2 large tare keys; clear and logically arranged keypad
- RS-232 serial interface port
- Four digital filter levels allow the balance to be optimally adapted to the particular ambient conditions at the place of installation.
- ISO/GLP-compliant logging of calibration, adjustment and weighing data using an optional data printer or a PC
- Built-in application programs: counting, weighing in %, net-total formulation; animal weighing; mass unit conversion by toggling between 2 units

Readability	0.001 g
Weighing capacity	320 g
Pan size	Ø 110* mm
Response time (avg.)	1.5 s
Housing (WxDxH)	213x342x153 mm

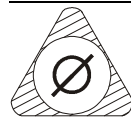

*) Triangular weighing pan shape. Ø = diameter of the inner circle; the shaded area is completely available for use



Data printer



SartoConnect

Key accessories	Order no.
Data printer with date, time and statistics functions	YDP03-OCE
Analytical draft shield chamber	YDS01CP
Dust cover	6960CP02
Data transfer software SartoConnect with RS-232C standard cable (1.5 m);	YSC01L
Calibration weight (1x200 g; F1)	YCW5238-00